



People.Health.Care.

Compression Therapy Toolbox

Products to support your patients for life.





People.Health.Care.



Achieve the best outcomes. Individualize your therapy.

L&R offers a full range of effective compression solutions, from the “gold standard” of bandaging to the Solaris Collection™ of garments and wraps for self-management.

When you choose L&R products, you can be confident your patients are in good hands with trusted brands like TributeNight™, ReadyWrap®, Tribute® Wrap, and Rosidal® K.

Our bandaging range offers the largest variety of products and sizes to best conform to the range of patient shape, size, and tissue needs you see every day.

The Solaris Collection™ — designed by lymphedema therapists through patient feedback — offers a range of products that supports lymphedema treatment from intensive phase all the way through self-management with daytime, nighttime, custom, and off-the-shelf solutions.

The Solaris Collection continues to expand to address your needs by listening to and understanding your clinical challenges and developing new products which leverage creative ideas about how to best meet your patients’ needs.

L&R is passionate about delivering innovative wound care and compression solutions for improved quality of life, as well as top-notch service to support you and your patients in their journey with lymphedema.



Solaris
Collection™

BY L&R USA INC.

Customer Service

Our Customer Service teams are known for their exceptional service, and can work with therapists, patients, and providers.

United States

🕒 8:00 AM–6:00 PM Central
Monday–Friday

☎ 1-855-892-4140

✉ inquiries@us.LRmed.com

Canada

🕒 8:00 AM–5:00 PM Eastern
Monday–Friday

☎ 1-888-836-6726

✉ inquiries@ca.LRmed.com

Custom Design Center

We offer therapists and fitters direct access to our team of designers, to lend their expertise in helping design custom garments for patients from coast to coast.

🕒 8:00 AM–6:00 PM Central, Monday–Friday

☎ 1-414-892-5158

✉ CustomDesignCenter@us.LRmed.com

📺 Video conference

TributeNight™

Tried and trusted for over 20 years.

TributeNight™ is a custom-made therapeutic foam garment for self-management during the evening, nighttime, and low activity.



Custom



Evening & Nighttime



Low Activity



Example garment styles





Why choose TributeNight™?

TributeNight™ was designed by lymphedema therapists to replace nightly multi-layer bandaging, offering patients an easy solution to maintain their therapy success on their own. Its comfortable, medium-stretch material provides containment to help the proven foam technology nestle into tissue to promote skin health¹ — key for lymphedema patient self-care.

TributeNight offers customized therapy, incorporating patient needs for style, fit, functionality, and donning preferences in every garment.




- 12-month warranty
- Easy garment care in washer and dryer

SURE-FIT™ GUARANTEE

Measuring for custom garments is a combination of art and science, and we want to ensure you are confident in the design and customized fit of each TributeNight.

Each TributeNight garment is covered, at no extra charge, for the life of the garment by our industry-leading SURE-FIT™ Guarantee. This complimentary guarantee is comprised of three industry-leading services:

- First Fit Warranty
- Reduction Alteration Warranty
- Materials & Manufacturing Warranty

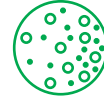

Proven foam technology







First-use chipped foam intended to create pressure differentials and promote a micro-massage-type effect to break up and soften hardened skin or fibrotic tissue to support skin health¹.



Powered by Coolcore®

Chemical-free Cooling Power rapidly wicks away sweat and moisture to regulate body temperature.

Building a custom TributeNight™ for your patient is easy.

-  **Determine garment style**
Coverage from head to toe and in-between.
-  **Select channeling type**
Designed based on MLD concepts.
-  **Choose profile height**
Offering foam technology variety benefiting patient comfort.
-  **Pick fabric color**
Patients can express themselves.
-  **Add modifications**
Nearly a dozen optional features available to enhance patient compliance.
-  **Add accessories**
Enhance the garment and their experience.



Helpful measuring videos are available on our website!



ReadyWrap®

Clinically-effective and easy-to-use.

ReadyWrap® is an off-the-shelf adjustable gradient compression wrap for transitioning out of therapy into self-management during the daytime, nighttime, and high or low activity.

  Off-the-Shelf or Custom*

 Day

  High or Low Activity



Available styles





Why choose ReadyWrap®?

ReadyWrap® features low-stretch straps that mimic short-stretch bandaging to provide clinically-effective compression therapy with high working and low resting pressure, all while supporting patient compliance with its unique alternating strapping system which can be donned with a single hand.


The modular design of the ReadyWrap family allows each garment to be worn on its own or in conjunction with other styles.

With ReadyWrap, you can be certain your patients will have an effective and comfortable garment that easily adjusts to their limb and donning abilities.




- Upper extremity styles: 20–40 mmHg
- Lower extremity styles: 30–50 mmHg

- 6-month warranty
- Hand wash, lay flat to dry

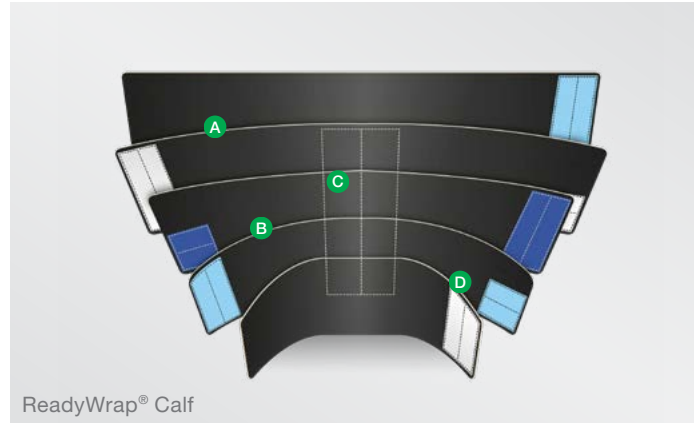

Low-stretch material

Mimics short-stretch bandaging — the gold standard of treatment.

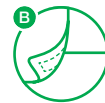

50% Overlapping straps[†]

Helps prevent gapping.²



Blocked support spine[‡]



Helps maintain garment form to support all-day wear.²



Color-coded straps

Provides easy step-by-step donning guidance.

Accessories

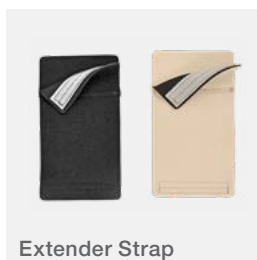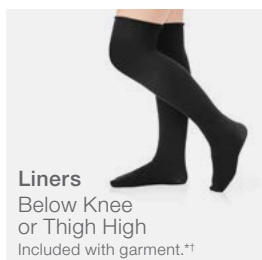



* Lower extremity styles only
† Except ReadyWrap Toe
‡ ReadyWrap Calf, Knee and Thigh styles only

Tribute® Wrap

The adjustability of ReadyWrap meets the legendary performance of TributeNight.

Tribute® Wrap is an off-the-shelf adjustable foam gradient compression garment for self-management during the evening, nighttime, and low activity.

 Off-the-Shelf

 Evening & Nighttime

 Low Activity



Available styles





Why choose Tribute® Wrap?

Tribute® Wrap is an innovative garment built with proven foam technology, which nestles into the tissue to promote skin health¹. Each garment is fully adjustable, while comfortably adjusting to your patient's unique shape and lifestyle. Its adjustability helps patients maintain gains made during therapy, while allowing the flexibility to support further reduction.

Simply pull it on, adjust the straps, and slide on the Sleep Sleeve for a comfortable compression experience.




Black exterior
Slate interior

- 6-month warranty
- Easy garment care in washer and dryer


Accessories



Knee to Thigh

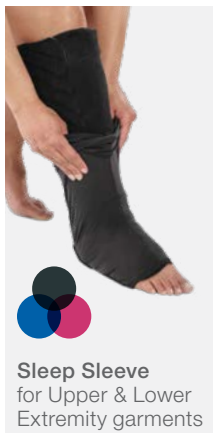


Sleep Sleeve for Upper & Lower Extremity garments

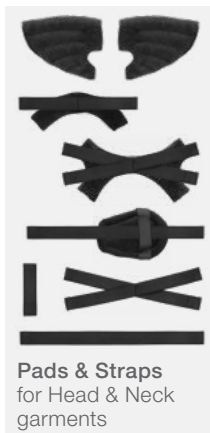

Pads & Straps for Head & Neck garments




Tribute® Wrap Wrist to Axilla

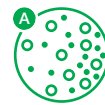

Proven foam technology

First-use chipped foam intended to create pressure differentials and promote a micro-massage-type effect to break up and soften hardened skin or fibrotic tissue to support skin health¹.

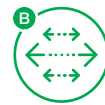

Adjustable strapping system

Easily conforms for the right, comfortable fit.



Easy-to-pull finger grip

Simple pull function to support self-management.




Powered by Coolcore®

Innovative interior fabric with chemical-free Cooling Power which rapidly wicks away sweat and moisture to regulate body temperature.

ExoStrong™

Affordable alternative to custom flat-knit.

ExoStrong™ is an off-the-shelf flat-knit gradient compression garment for self-management during the daytime or high activity.

 Off-the-Shelf

 Daytime

 High activity



Available styles






Why choose ExoStrong™?

For lymphedema patients, flat-knit garments are one of the primary self-care options. Where circular-knit material can migrate into skin folds, flat-knit helps shape the limb and provides the high containment needed to maintain reductions made in therapy. Additionally, ExoStrong™'s soft, breathable design helps support patient compliance.

ExoStrong is a fraction of the cost of custom flat-knit garments, ensuring patients don't have to choose between an effective maintenance garment and one they can afford.

up to
60%
savings
compared to
custom flat knit

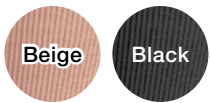

Compression	±† 20–30 mmHg
Top style	KNIT TOP [*] ▨ SILICONE TOP [†]
Toe style	○ OPEN TOE
Digit style	☞ QUARTER HALF FULL

^{*} Knit top is available with Arm and Below Knee
[†] Silicone top is available with Arm, Below Knee, and Thigh High

- 6-month warranty
- Easy garment care in washer and dryer



Made in the USA

High-quality flat knit manufactured at our headquarters in Milwaukee, WI.



Supports mobility

High containment garment in a comfortable and functional design.



Quick and easy fitting

Only 2–5 measurements needed to select size for each style.



Ships quickly



Orders placed before 2 PM Central / 3 PM Eastern ship the same day.



SwellSpots®

Small pads that offer big relief!

SwellSpots® are reusable therapeutic foam pads to incorporate into multi-layer bandaging or garments to provide focal compression to help soften troublesome areas of edema and fibrosis.

  Off-the-Shelf or Custom

  Day or Night

  High or Low Activity



Available off-the-shelf styles





“My patients find SwellSpots not only make it really easy to break up tough fibrosis and scar tissue, but also help alleviate any pain — making it easier to be compliant with their garments.”

Maureen McBeth, PT, CLT-LANA




Why choose SwellSpots®?

Using these versatile pads will save time spent cutting your own shapes or creating chip foam bags to allow more time for manual lymphatic drainage (MLD), bandaging, and patient education.

Select a pad style and size from our wide range of options, place either side of the pad directly against the skin then apply short stretch bandages over the top, or tuck into a garment. From there, let the pad do the rest.

Whether it's breaking up fibrotic tissue, filling concavities, or padding bony prominences, these small SwellSpots® pads offer you and your patients a big relief!

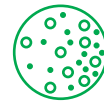

Proven foam technology

First-use chipped foam intended to create pressure differentials and promote a micro-massage-type effect to break up and soften hardened skin or fibrotic tissue to support skin health¹.



Versatile shapes & sizes

Wide range of off-the-shelf and custom solutions to manage troublesome areas no matter the location.



Channeling

Designed based on MLD principles of directing fluid along lymphatic pathways.



Reusable pads

Durable construction ensures pads retain their shape between washes.




Slate

- 6-month warranty for off-the-shelf pad styles
1-year warranty for custom pads
- Easy care in washer and dryer

Caresia™

Simplifies and enhances your bandaging protocol.

Caresia™ is an off-the-shelf foam bandage liner designed to be used during the intensive phase of therapy to simplify the multi-layer bandaging process during the daytime, nighttime, and high or low activity.

-  Off-the-Shelf
-  Day or Night
-  High or Low Activity



Available styles



Gauntlet



Wrist to Axilla



MCP to Axilla



Foot



Thigh



Below Knee



Glove



Before CDT



After CDT

Why choose Caresia™?

Caresia™ replaces the underlayers of standard multi-layer bandaging protocol, including the skin protection layer (e.g., stockinette), padding, and focal compression/foam.

With its high stretch fabric, Caresia features great conformity, even around hard-to-manage lobules or with unique presentations.

Simply slide it on and bandage over the top with short-stretch bandages like Rosidal® K to achieve the desired compression profile.



Saves time

Eliminates the underlayers of bandaging. Slide on and apply short-stretch bandaging over top.



Proven foam technology

First-use chipped foam intended to create pressure differentials and promote a micro-massage-type effect to break up and soften hardened skin or fibrotic tissue to support skin health¹.



Quick and easy fitting

Only 2–5 measurements needed to select size for each style.



High-stretch fabric

Offers great range of motion compared to multi-layer bandages and conforms well to body contours.

Accessories



Rosidal® K
Short-stretch bandage






- 6-month warranty
- Easy garment care in washer and dryer

ExoSoft™

Known for its soft, comfortable, and high containment hand garments.

ExoSoft™ is an off-the-shelf circular-knit gradient compression garment designed for self-management during the day or high activity.

-  Off-the-Shelf
-  Day
-  High Activity



- 6-month warranty
- Easy garment care in washer and dryer

Compression	15–20 mmHg*	20–30 mmHg
Top style	KNIT TOP	▨ SILICONE TOP†
Toe style	○ OPEN TOE	● CLOSED TOE
Digit style	☞ QUARTER	HALF FULL

* 15–20 mmHg is available with Below Knee and Thigh High
 † Silicone top is available with Arm Sleeve and Thigh High

Available styles



ExoCustom™



Custom made in the USA.

ExoCustom™ is a custom flat-knit compression garment for self-management during the day or high activity.

 Custom

 Day



 High Activity

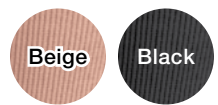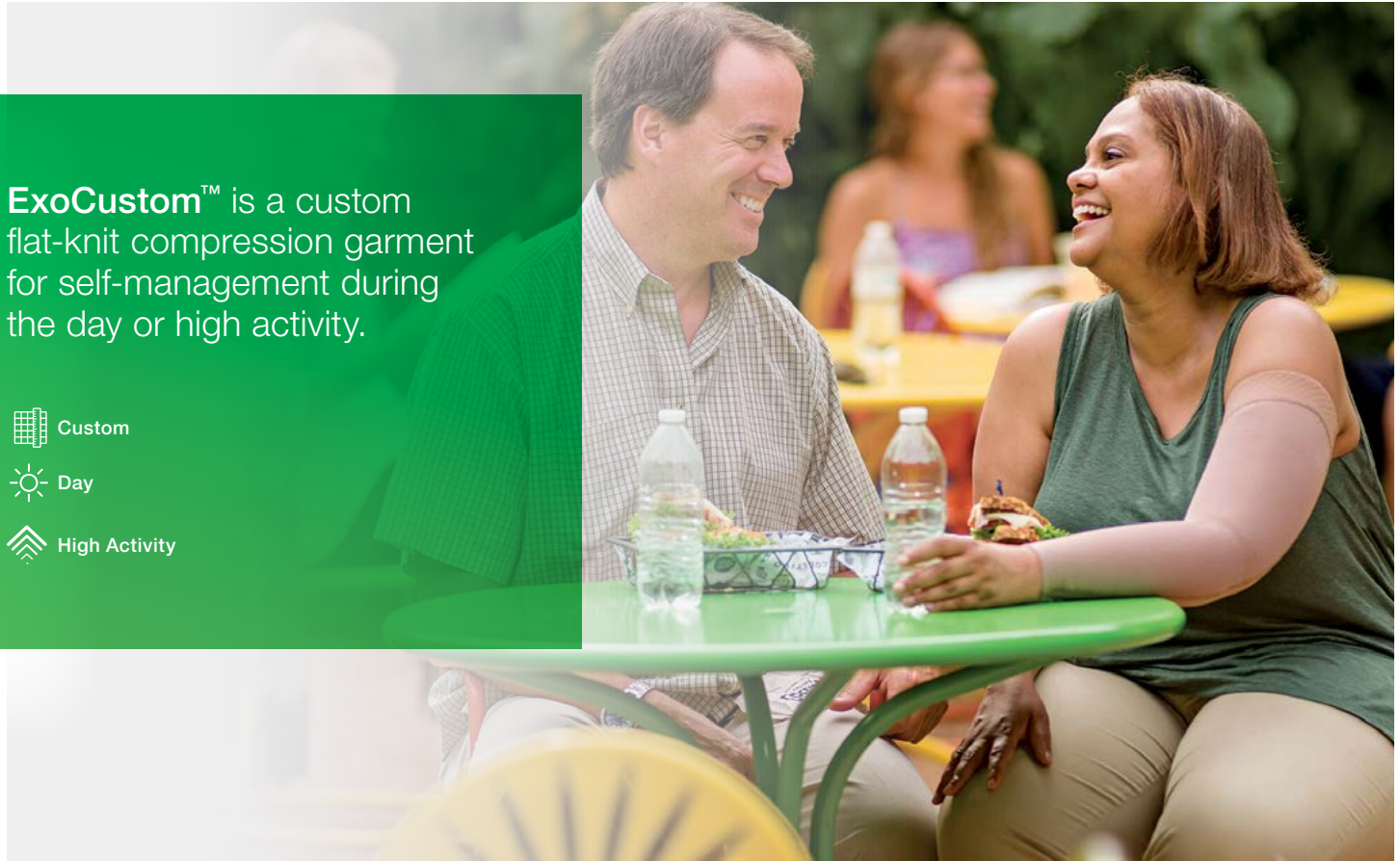

- 6-month warranty
- Easy garment care in washer and dryer

Compression	↕ 18–21 mmHg	23–32 mmHg	34–46 mmHg	
Silicone	WIDTH	DIAMETER	COVERAGE	PLACEMENT LOCATION
Pocket	ANY PLACEMENT LOCATION			
Zipper	ANY PLACEMENT LOCATION			

Available styles

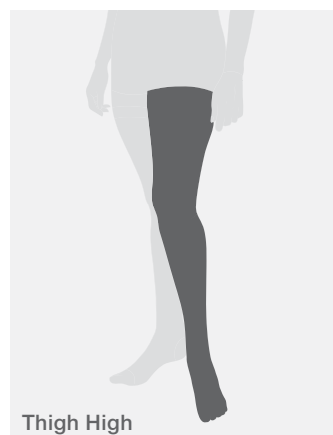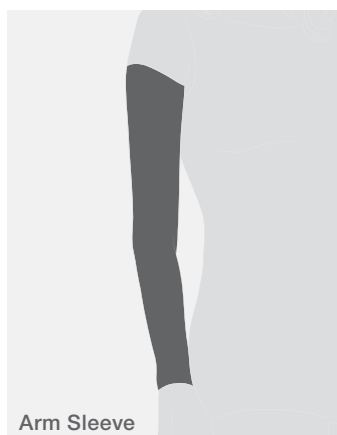

Multi-layer bandaging products

The widest range of bandage types and sizes available

STEP 1

Skin Protection



tg® tubular stockinette



tg® cotton stockinette



tg® soft stockinette

STEP 2

Fingers/Toes



Transelast® Classic conforming bandage



Mollelast® conforming bandage



Raucolast® conforming bandage

STEP 3

Padding



Cellona® synthetic padding



Rosidal® soft foam padding



Caresia™ foam bandage liner



lable domestically.

STEP 4

Focal Compression



Komprex®
medical-grade foam pad, roll, and sheet



Komprex® II
channel foam sheet



SwellSpots™
edema and fibrosis pad



STEP 5

Bandaging



Rosidal® K
short-stretch bandage



Idealbinde®
wide short-stretch bandage



Lenkelast®
medium-stretch bandage



Durelast®
very-short-stretch bandage

STEP 6

Tape



Silkafix®
adhesive tape



People.Health.Care.



References

- 1. Chohan, A., Haworth, L., Sumner, S., Olivier, M., Birdsall, D., and Whitaker, J. (2019). Examination of the effects of a new compression garment on skin tissue oxygenation in healthy volunteers. *Journal of Wound Care*. 28. 1.
- 2. Ehmann, S., and Bock, K. (Jan. 2018), "Case report to demonstrate the need for selection criteria for optimal adjustable Velcro wrap prescription", *The Journal of Advanced Wound Care, North American Supplement*, vol. 27, no. 1

Coolcore® is a registered trademark of Coolcore, LLC. VELCRO® is a registered trademark of Velcro BVBA. CertiPUR-US® is a certification program administered by a not-for-profit organization. Although the CertiPUR-US® foam certification program is limited to and exclusively for flexible polyurethane foams used in bedding, pillows and upholstered furniture, our foams pass the same rigorous standards set by the CertiPUR-US® program.